



4 pillars of design

1 USER COMFORT

Utilising Dreamline's Bubble Contour design for comfortable immersion with breathability channels. A central channel is included for spinal relief. The foam pad is overlayed with a pressure redistributing viscoelastic "Sunmate" foam for extra long-term comfort.

FLO2™ engineered 4-way ultra-stretch fabrics with moisture wicking properties enable the user to fully utilise the comfort of the pad.


The pad & cover lock to the shell securely for trouble free riding.

The cover is removable & machine washable.

2 CONTOUR PROFILE

Specifically designed with an open curvature shape with a more closed lateral angle to encourage mid-line sitting posture.

Slim 1" wide laterals maximise free arm movement & ensure a compact user-to-wheelchair fit.

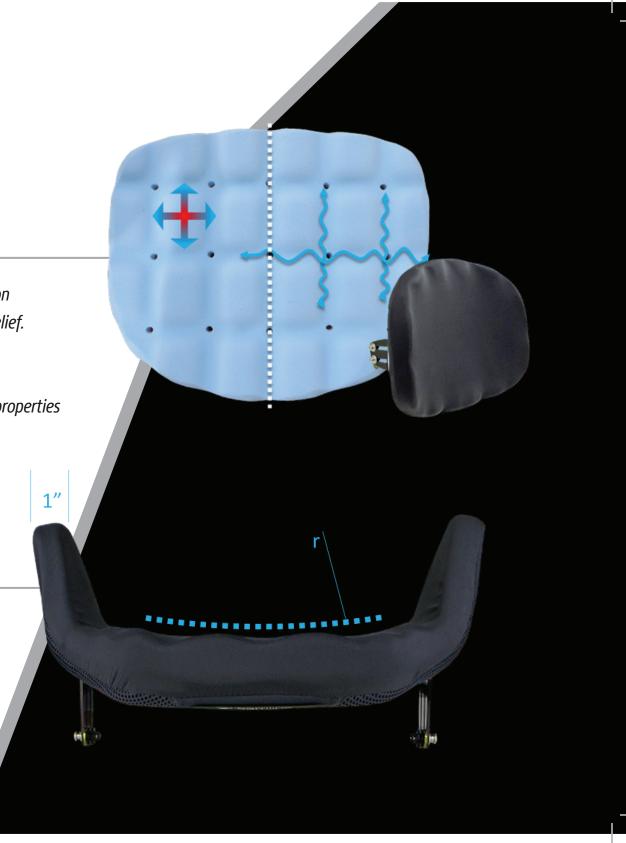



3 POSITIONING

Ignite backrests have a series of positioning pads supplied in the standard kit including a Lumbar pad & a pair of wide wedges & narrow wedges.

The added ability to include Propoline's BAC Backrest Advanced

The added ability to include Dreamline's BAC Backrest Advanced
Contouring system to any Ignite backrest ensures users with
asymmetrical postures will be well supported.
Truly reversible with an offset lateral design ensuring the
lateral height best suits the user.

4 ADJUSTABILITY

The proven U2 quick-release hardware is effortless & quick to operate. The U2
Compact clamps attach to tubing diameters from 19mm (3/4") to 25mm (1").

Enormous width adjustability (4") ensures the right fitting backrest is selected for the user & ensures it will fit the wheelchair frame.

Up to 60deg. of angle adjustment (+/- 30deg.) & 5" of depth adjustment takes the backrests ability beyond the need.



eding(

easy selection

2 types


CARBON

ALUMINIUM

3 contour depths





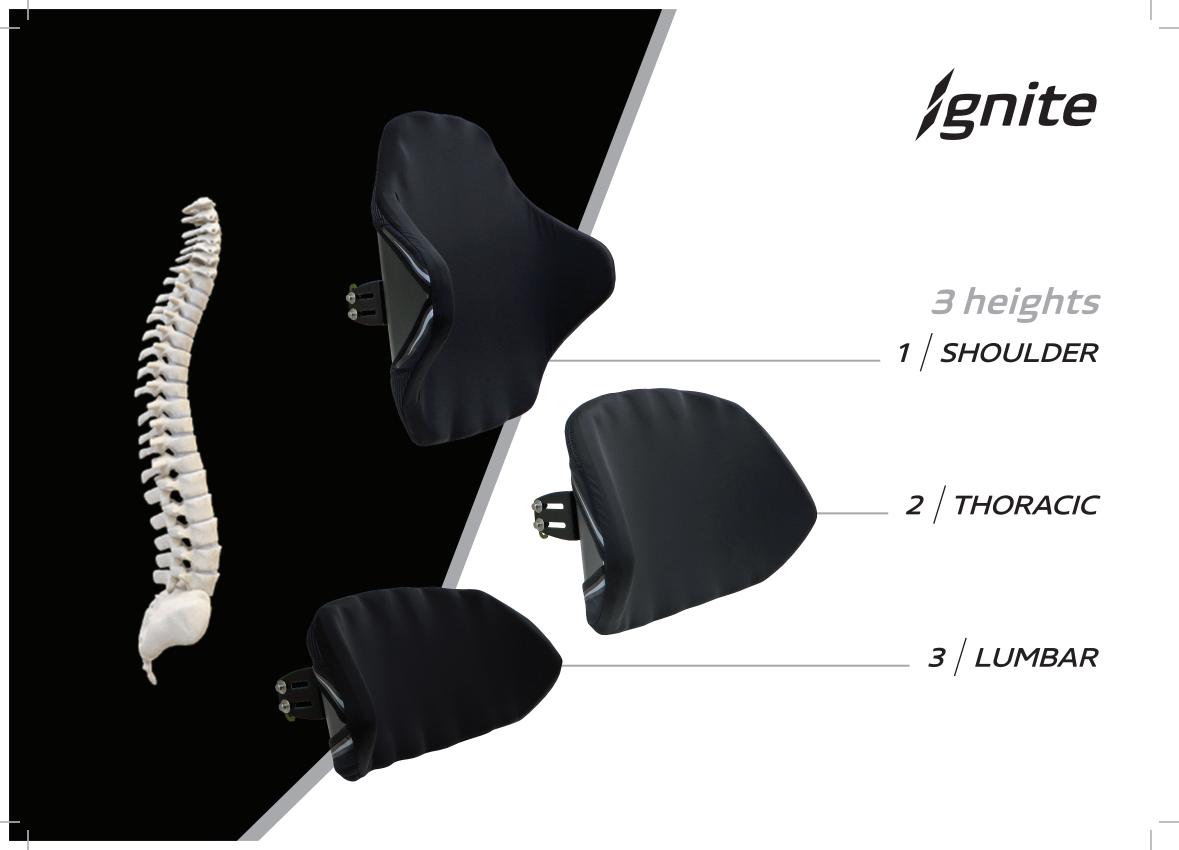

/gnite pro carbon series




Dreamline Ignite Pro backrests have a lightweight solid weave carbon fibre shell. This makes for a much stronger shell than hollow core or honeycomb designs.


The weave pattern has been carefully designed to enable a limited flex for comfort with maximum strength for every-day use.

adjustable aluminium laterals



Dreamline **Ignite Pro** Backrests with Mid or Deep Contour options incorporate robust aluminium laterals which attach to the shell through the same fixings as the backrest mounting brackets for added strength.

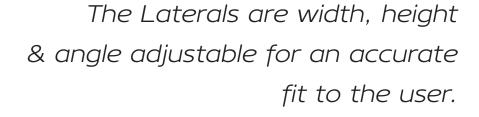

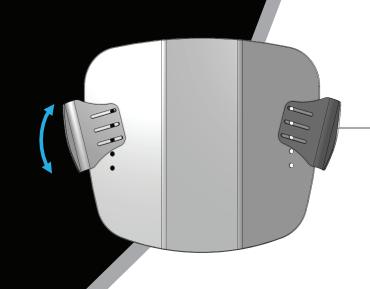

active / mid / deep contours

empowering performance

with the right backrest, postural efficiency & sitting ability can increase, helping to improve overall performance

/gnite



accessories + upgrades



BAC System for asymmetrical postures



Compact L-brackets



Extra Depth Kit





4-point Upgrade (Tall Backs)



Lateral Stability Wedge



Lumbar Pad



External Swing-away Laterals (Active Contour backs)



Spare Outer Cover




Incontinence Cover




gnite

COLOUR HIGHLIGHTS

"Spacer" fabric Panels on the border are swapped out with colour panels.

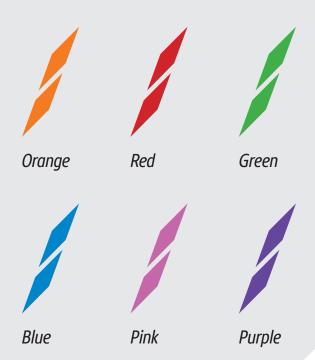

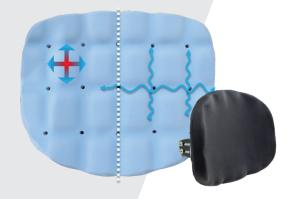

FLO2™ ENGINEERED FABRIC

FLO2[™] is an advanced high-performance fabric designed for medical and high energy sporting users.











Extreme
4-Way Stretch

Air Exchange

Moisture Wicking

Durable

Fire Retardant



tested for safety


The **Dreamline Ignite Backrests U2**

compact hardware is fully crash tested to meet the standards of WC20 & ISO7176-19 for forward facing transportation with an occupant in a motor vehicle.

When occupied in a motor vehicle, Dreamline backrests must only be used on wheelchair frames which provide four anchorage points & pelvic belt anchor points which comply with the requirements of RESNA WC-4:2012, Section 19.



Rolapal Ltd 03 4135 2558 sales@rolapal.com.au www.rolapal.co.nz

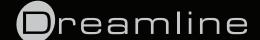

gnite

R•LAPAL